

Smart-IP Series

Anti-V5-tag mAb-Magnetic Beads

CODE No.	M167-11
CLONALITY	Monoclonal
CLONE	1H6
ISOTYPE	Mouse IgG2a κ
QUANTITY	20 tests (Slurry: 1 mL)
SOURCE	Purified IgG from hybridoma supernatant
IMMUNOGEN	Carrier protein conjugated synthetic peptide, GKPIPPLLGLDST (V5-tag)
FORMULATION	10 mg magnetic beads in 1 mL PBS/0.1% BSA/0.09% NaN ₃
	*Azide may react with copper or lead in plumbing system to form explosive metal azides. Therefore, always flush plenty of water when disposing materials containing azide into drain.
STORAGE	This gel slurry is stable for one year from the date of purchase when stored at 4°C.

APPLICATION-CONFIRMED

Immunoprecipitation 50 μ L of beads slurry/sample

*The purification capacity of Anti-V5-tag mAb-Magnetic Beads varies depending upon the characteristics of a V5-tagged protein.
For example, 50 μ L of beads slurry bounds 3.2 μ g of a V5-tagged protein (32 kDa).

APPLICATION-REPORTED


Co-Immunoprecipitation 50 μ L of beads slurry/sample

For more information, please visit our website <https://ruo.mbl.co.jp/>.

The descriptions of the following protocols are examples. Each user should determine the appropriate condition.

Immunoprecipitation

- 1) Wash 5×10^6 cells 3 times with PBS and suspend with 1 mL of cold Lysis buffer [50 mM Tris-HCl (pH 7.5), 150 mM NaCl, 1% NP-40] containing appropriate protease inhibitors.
- 2) Centrifuge the tube at $12,000 \times g$ for 5 min. at 4°C and transfer the supernatant to another tube.
- 3) Add magnetic beads as suggested in the **APPLICATION** into 400 μL of the supernatant prepared in step 2). Mix well and incubate with gentle agitation for 1 hr. at 4°C .
- 4) Place the tube on the magnetic rack (MBL; code no. 3190) for a few seconds.
- 5) Remove the supernatant.
- 6) Wash the beads 4 times with 1 mL of cold Wash buffer [50 mM Tris-HCl (pH 7.5), 150 mM NaCl, 0.05% NP-40] (place the tube on the magnetic rack for a few seconds).
- 7) Resuspend the magnetic beads in 20 μL of Laemmli's sample buffer, boil for 2 min., and place the tube on the magnetic rack for a few seconds.
- 8) Load 20 μL of the sample per lane in a 1-mm-thick SDS-polyacrylamide gel and carry out electrophoresis.
- 9) Visualize the protein bands by CBB staining.

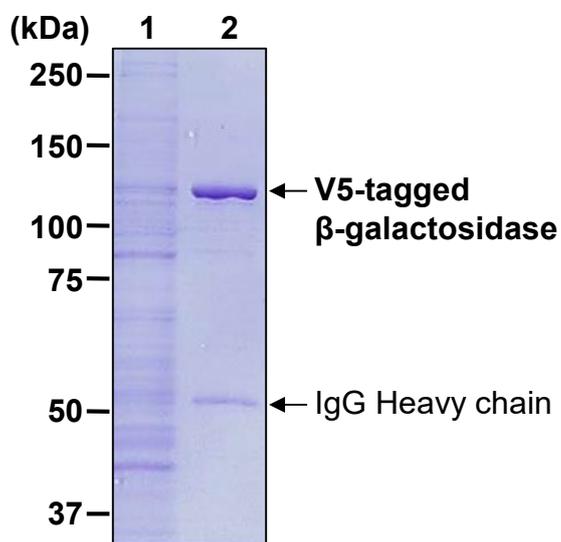

Immunoprecipitation of V5-tagged protein

Sample: V5-tagged β -galactosidase/HEK293T whole cell lysate

Lane 1: Input (5 μL /lane)

Lane 2: Post-IP beads of Anti-V5-tag mAb (MBL; code no. M167-11)