

 **My select** sampler set

# Anti-NRF2 pAb

**CODE No.** PM069MS

**CLONALITY** Polyclonal  
**ISOTYPE** Rabbit Ig, affinity purified  
**QUANTITY** 20 µL

**SOURCE** Purified Ig from rabbit serum  
**IMMUNOGEN** Human NRF2, full length (recombinant)  
**FORMULATION** PBS containing 50% Glycerol (pH 7.2). No preservative is contained.  
**STORAGE** This antibody solution is stable for one year from the date of purchase when stored at -20°C.

## APPLICATIONS-CONFIRMED

Western blotting 1:1,000  
Immunoprecipitation 5 µL/300 µL of cell extract from 3 x 10<sup>6</sup> cells  
Immunocytochemistry 1:1,000  
Immunohistochemistry 1:1,000 (paraffin section)  
Heat treatment for paraffin embedded section: microwave oven, for 20 min. in 10 mM citrate buffer (pH 6.0)

## SPECIES CROSS REACTIVITY on WB

Species	Human	Mouse	Rat	Hamster
Cells	Transfectant, HeLa, A549	NIH/3T3, WEHI-3B, L5178Y, C2C12	PC12	CHO
Reactivity	+	+ (weak)	+ (weak)	+ (weak)

**Entrez Gene ID** 4780 (Human)

**REFERENCES**


- 1) Taguchi, K., *et al.*, *Genes Cell* **16**, 123-140 (2011)
- 2) Komatsu, M., *et al.*, *Nat. Cell Biol.* **12**, 213-223 (2010)
- 3) Nguyen, T., *et al.*, *J. Biol. Chem.* **284**, 13291-13295 (2009)

For more information, please visit our web site <https://ruo.mbl.co.jp/>.

**SDS-PAGE & Western blotting**

- 1) Wash  $1 \times 10^7$  cells 3 times with PBS and suspends them in 1 mL of Laemmli's sample buffer, then sonicate briefly (up to 20 sec.).
- 2) Boil the samples for 3 min. and centrifuge. Load 10  $\mu$ L of the sample per lane in a 1-mm-thick SDS-polyacrylamide gel (7.5% acrylamide) for electrophoresis.
- 3) Blot the protein to a polyvinylidene difluoride (PVDF) membrane at 1 mA/cm<sup>2</sup> for 1 hr. in a semi-dry transfer system (Transfer Buffer: 25 mM Tris, 190 mM glycine, 20% methanol). See the manufacturer's manual for precise transfer procedure.
- 4) To reduce nonspecific binding, soak the membrane in 10% skimmed milk (in PBS, pH 7.2) overnight at 4°C.
- 5) Wash the membrane with PBS-T (0.05% Tween-20 in PBS) (5 min. x 3).
- 6) Incubate the membrane with primary antibody diluted with 1% skimmed milk (in PBS, pH 7.2) as suggested in the **APPLICATIONS** for 1 hr. at room temperature. (The concentration of antibody will depend on the conditions.)
- 7) Wash the membrane with PBS-T (5 min. x 3).
- 8) Incubate the membrane with the 1:10,000 Anti-IgG (Rabbit)-HRP (MBL, code no. 458) diluted with 1% skimmed milk (in PBS, pH 7.2) for 1 hr. at room temperature.
- 9) Wash the membrane with PBS-T (5 min. x 3).
- 10) Wipe excess buffer on the membrane, and then incubate it with appropriate chemiluminescence reagent for 1 min. Remove extra reagent from the membrane by dabbing with paper towel, and seal it in plastic wrap.
- 11) Expose to an X-ray film in a dark room for 1 min. Develop the film as usual. The condition for exposure and development may vary.

(Positive controls for Western blotting; Transfectant, HeLa, A549)

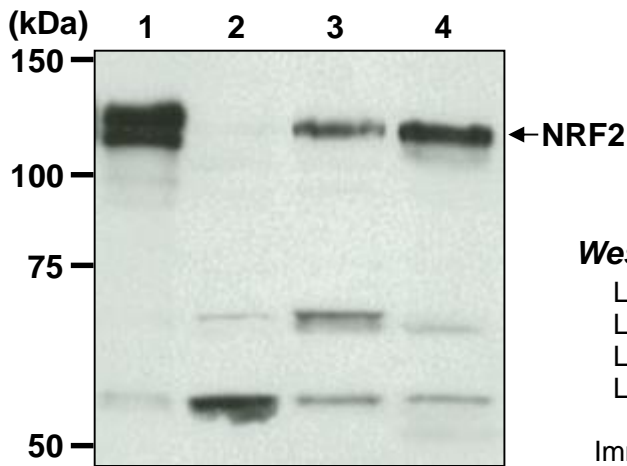

***Western blot analysis of NRF2***


Lane 1: NRF2/293T  
Lane 2: 293T  
Lane 3: HeLa  
Lane 4: A549

Immunoblotted with Anti-NRF2 pAb (MBL, code no. PM069)

### Immunoprecipitation

- 1) Wash  $1 \times 10^7$  cells twice with PBS and resuspend them with 1 mL of ice-cold Lysis buffer [50 mM Tris-HCl (pH 7.5), 150 mM NaCl, 0.05% NP-40] containing appropriate protease inhibitors, then sonicate briefly (up to 20 sec.).
- 2) Centrifuge the tube at 12,000 x g for 5 min. at 4°C and transfer the supernatant to another tube.
- 3) Mix 20  $\mu$ L of 50% protein A agarose beads slurry resuspended in 300  $\mu$ L of IP buffer [10 mM Tris-HCl (pH 8.0), 500 mM NaCl, 0.1% NP-40] with primary antibody as suggested in the **APPLICATIONS**. Incubate with gently agitation for 1 hr. at room temperature.
- 4) Wash the beads 3 times with 1 mL of IP buffer.
- 5) Add 300  $\mu$ L of cell lysate (prepared sample from step 2)), then incubate with gentle agitation for 1 hr. at room temperature.
- 6) Wash the beads 5 times with 1 mL of Lysis buffer.
- 7) Resuspend the beads in 20  $\mu$ L of Laemmli's sample buffer, boil for 3 min. and centrifuge.
- 8) Load 10  $\mu$ L of the sample per lane in a 1-mm-thick SDS-polyacrylamide gel (7.5% acrylamide) for electrophoresis.
- 9) Blot the protein to a polyvinylidene difluoride (PVDF) membrane at 1 mA/cm<sup>2</sup> for 1 hr. in a semi-dry transfer system (Transfer Buffer: 25 mM Tris, 190 mM glycine, 20% methanol). See the manufacturer's manual for precise transfer procedure.
- 10) To reduce nonspecific binding, soak the membrane in 10% skimmed milk (in PBS, pH 7.2) overnight at 4°C.
- 11) Wash the membrane with PBS-T (0.05% Tween-20 in PBS) (5 min. x 3).
- 12) Incubate the membrane with 1:1,000 Anti-NRF2 pAb (MBL, code no. PM069) diluted with 1% skimmed milk (in PBS, pH 7.2) for 1 hr. at room temperature. (The concentration of antibody will depend on the conditions.)
- 13) Wash the membrane with PBS-T (5 min. x 3).
- 14) Incubate the membrane with the 1:10,000 Anti-IgG (Rabbit)-HRP (MBL, code no. 458) diluted with 1% skimmed milk (in PBS, pH 7.2) for 1 hr. at room temperature.
- 15) Wash the membrane with PBS-T (5 min. x 3).
- 16) Wipe excess buffer on the membrane, and then incubate it with appropriate chemiluminescence reagent for 1 min.
- 17) Remove extra reagent from the membrane by dabbing with paper towel, and seal it in plastic wrap.
- 18) Expose to an X-ray film in a dark room for 1 min. Develop the film as usual. The condition for exposure and development may vary.

(Positive control for Immunoprecipitation; HeLa)

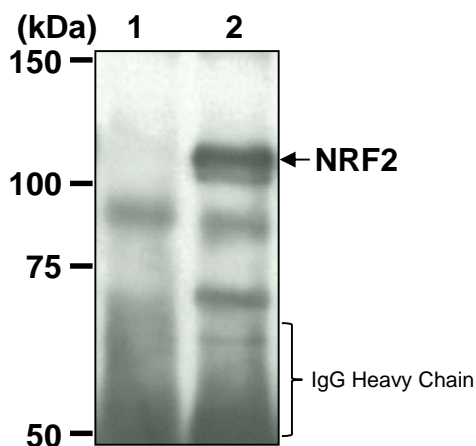

#### **Immunoprecipitation of NRF2 from HeLa**

Lane 1: IP with normal rabbit IgG (MBL, code no. PM035)


Lane 2: IP with Anti-NRF2 pAb (MBL, code no. PM069)

Immunoblotted with Anti-NRF2 pAb (MBL, code no. PM069)

### **Immunocytochemistry**

- 1) Spread the cells on a glass slide, and then incubate in a CO<sub>2</sub> incubator overnight.
- 2) Remove the culture supernatant by careful aspiration.
- 3) Wash the slide twice with PBS.
- 4) Fix the cells with 4% paraformaldehyde/PBS for 10 min. at room temperature (20~25°C).
- 5) Wash the slide twice with PBS.
- 6) Permeabilize the cells with 200 µL of 0.2% Triton X-100/PBS for 10 min. at room temperature.
- 7) Wash the slide twice with PBS.
- 8) Add 200 µL of the primary antibody diluted with 2% fetal calf serum (FCS)/PBS as suggested in the **APPLICATIONS** onto the cells and incubate for 1 hr. at room temperature. (Optimization of antibody concentration or incubation condition is recommended if necessary.)
- 9) Wash the slide twice with PBS.
- 10) Add 200 µL of 1:500 anti-IgG (Rabbit)-Alexa Fluor<sup>®</sup> 488 (Thermo Fisher Scientific, code no. A-11008) diluted with 2% fetal calf serum (FCS)/PBS onto the cells. Incubate for 30 min. at room temperature. Keep out light by aluminum foil.
- 11) Wash the slide twice with PBS.
- 12) Promptly add mounting medium onto the slide, then put a cover slip on it.

(Positive control for Immunocytochemistry; HeLa)

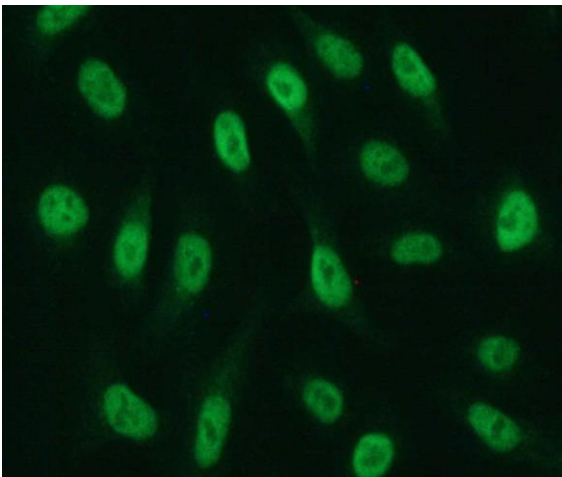


### ***Immunocytochemical detection of NRF2 in HeLa***

Green: Anti-NRF2 pAb (MBL, code no. PM069)

**Immunohistochemistry for formalin fixed paraffin-embedded section**

- 1) Deparaffinize the sections with Xylene 3 times for 3 min. each.
- 2) Wash the slides with Ethanol 3 times for 3 min. each.
- 3) Wash the slides with PBS 3 times for 3 min. each.
- 4) Remove the slides from PBS and heat-treated twice with 10 mM citrate buffer (pH 6.0) for 10 min. each using microwave.
- 5) Let the slides cool down at room temperature in the citrate buffer.
- 6) Wash the slides with running water for 5 min., then wash with PBS for 5 min.
- 7) Remove the slides from PBS and inactivate endogenous peroxidase with 3% H<sub>2</sub>O<sub>2</sub> in PBS for 10 min.
- 8) Wash the slides with running water for 5 min., then wash with PBS for 5 min.
- 9) Remove the slides from PBS, wipe gently around each section and cover tissues with blocking buffer [20 mM HEPES (pH 7.4), 1% BSA, 135 mM NaCl] for 5 min. at room temperature to block non-specific staining. Do not wash.
- 10) Tip off the blocking buffer, wipe gently around each section and cover tissues with primary antibody diluted with blocking buffer as suggested in the **APPLICATIONS**. (The concentration of antibody will depend on the conditions.) Incubate the sections for 1 hr. at room temperature.
- 11) Wash the slides 3 times in PBS for 5 min. each.
- 12) Wipe gently around each section and cover tissues with Histostar (Ms + Rb) (MBL, code no. 8460). Incubate for 30 min. at room temperature.
- 13) Wash the slides 3 times in PBS for 5 min. each.
- 14) Visualize by reacting for 10 min. with Histostar™ DAB Substrate Solution (MBL, code no. 8469). \*DAB is a suspect carcinogen and must be handled with care. Always wear gloves.
- 15) Wash the slides in water for 5 min.
- 16) Counter stain in hematoxylin for 1 min., wash the slides 3 times in water for 5 min. each, and then immerse the slides in PBS for 5 min.
- 17) Dehydrate by immersing in Ethanol 3 times for 3 min. each, followed by immersing in Xylene 3 times for 3 min. each. Now ready for mounting.

(Positive controls for Immunohistochemistry; human lung carcinoma, human colon carcinoma)

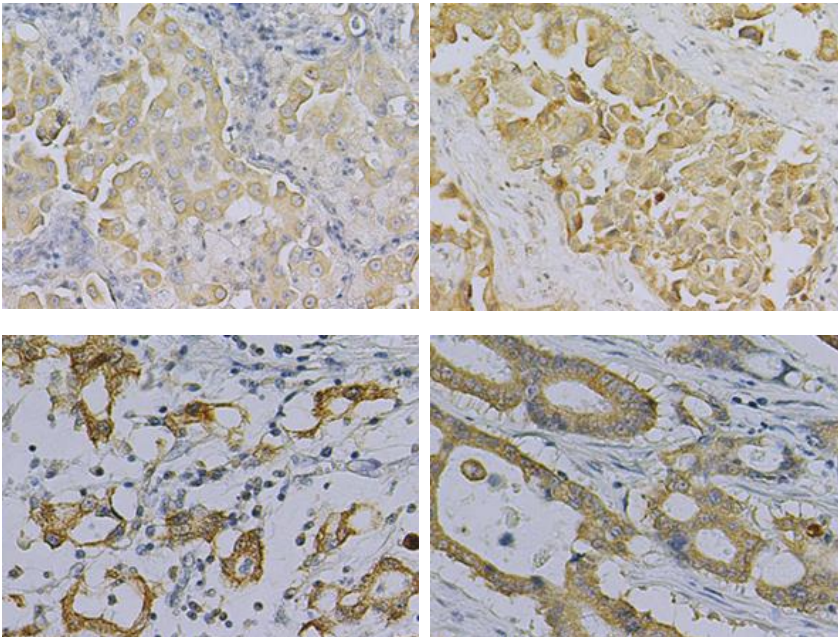

***Immunohistochemical detection of NRF2 in human cancer tissue***

Upper: Lung carcinoma (different fields)  
Lower: Colon carcinoma (different fields)